The BabyChecker™ Scale
La báscula BabyChecker™
Le pèse-bébé BabyChecker^{mc}


Assembly, Use and Service Manual

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The BabyChecker™ Scale

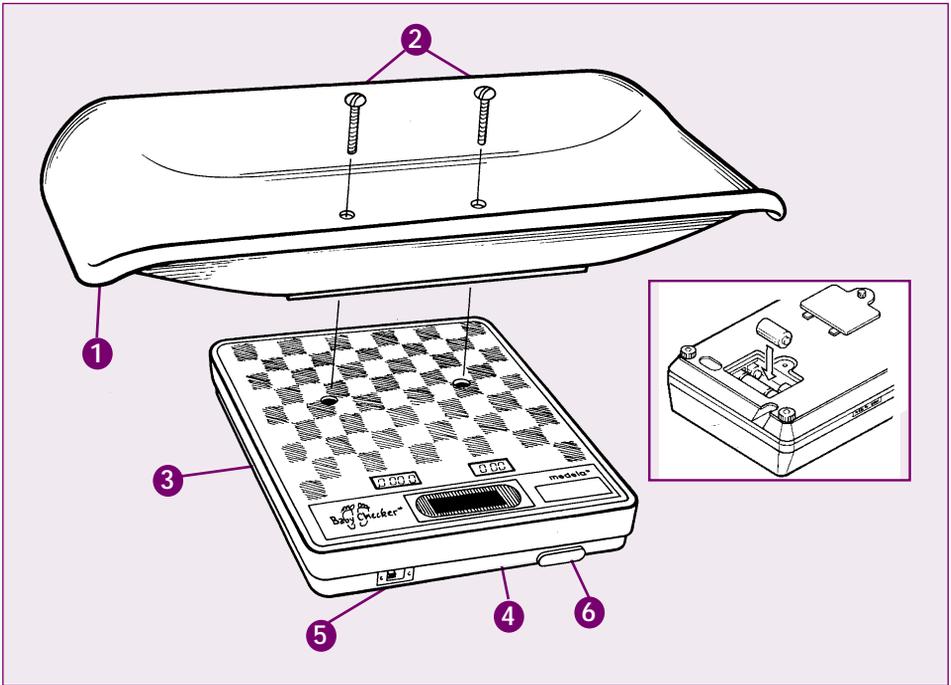Components

1. Baby Basket
2. 2 Fixing Screws
3. Scale Housing
4. Display Panel
5. lb•oz / off / kg Switch
6. Start/Zero Button
7. Battery Compartment Cover
8. Battery Compartment
9. 6 Alkaline AA Batteries

SETTING UP THE SCALE

■ Unpacking

Remove Scale Housing, Baby Basket and batteries from their packing materials.

■ Installing Batteries

The scale is powered by 6 AA alkaline batteries.

1. Remove Battery Compartment Cover, located on the bottom of the Scale Housing.
2. Install 6 AA batteries in the positions indicated on the inside bottom of the Battery Compartment.
3. Replace Battery Compartment Cover.


■ Assembly

The scale is packaged in two easy to assembly pieces.

If using the scale to weigh infants, the scale must be assembled before use. Follow these directions for assembly.

1. Remove Fixing Screws from Scale Housing.
2. Place Baby Basket on top of the Scale Housing. Align openings on Basket with openings on Scale Housing. Affix with 2 Fixing Screws.

If using the scale to weigh toddlers, the Baby Basket may be put aside and the scale housing used alone.

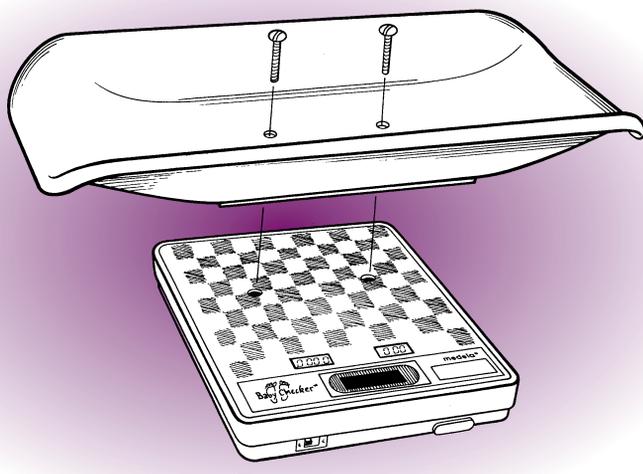■ Before Use

1. Choose a location for the scale. When lifting or moving the scale, be sure to lift the entire scale from the bottom of the Scale Housing, not the Baby Basket. To get an accurate weight, place the scale away from air conditioning and heat vents and direct sunlight. Place the scale on a stable, flat, level surface. **NEITHER THE BABY BASKET OR SCALE HOUSING SHOULD HANG OVER THE EDGE OF ANY SURFACE.**
2. Check the calibration of the scale following the instructions provided in the Service Section of this manual.

ZERO FUNCTION

The BabyChecker also has the feature of allowing a blanket, towel or other type of padding or protection to be placed on the scale before weighing and the weight of this piece being subtracted from the weight calculation.


To use this feature, place the towel, etc. on the Baby Basket or scale housing and then press the Start/Zero switch.

WEIGHING INFANTS

**WARNING: DO NOT LEAVE INFANTS UNATTENDED.
A FALL FROM THE SCALE COULD CAUSE SERIOUS INJURY.**

Before weighing infants, make sure the Baby Basket has been attached securely to the Scale Housing. The Baby Basket has been designed to accept weights up to 20 pounds or 10 kilograms. For heavier babies, read the section which follows on “Weighing Toddlers”.

The figures below indicate the method of the display.

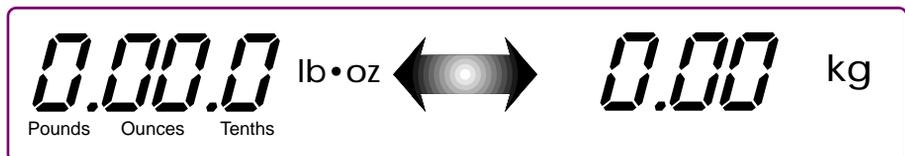
1. Turn the power on by sliding the lb•oz / off / kg switch to either pounds or kilograms, as you prefer.
2. Push the Start/Zero button.
3. If a blanket or other basket cover is to be used, place it on the scale. When the display stops fluctuating, press the Start/Zero button again.
4. When the display shows zero, place infant on the Baby Basket within 60 seconds.
5. When the display begins flashing, weight calculation is complete.
6. Remove infant from scale and then record weight. The weight will remain on the display for two minutes.
7. Although the BabyChecker will automatically turn off after about two minutes, slide the lb•oz / off / kg switch to OFF. Then if the Start/Zero button is accidentally pressed, the batteries will not be used.

WEIGHING TODDLERS

The unique design of the BabyChecker allows the Baby Basket to be removed and used similarly to a regular bathroom scale. When weighing babies between 20 – 40 pounds or 10 – 20 kilograms, be sure the Baby Basket is not in place on the scale.

The figures below indicate the method of the display.

1. Turn the power on by sliding the lb•oz / off / kg switch to either pounds or kilograms, as you prefer.

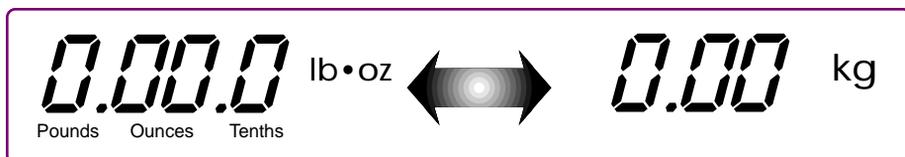2. Push the Start/Zero button.
3. When the display shows zero, place toddler on the Scale Housing within 60 seconds.
4. When the display begins flashing, weight calculation is complete.
5. If OL appears in the display, the total weight exceeds 20 kg or 40lb. An alternate scale will be necessary to determine actual weight.

6. Remove toddler from the scale and then record weight. The weight will remain on the display for two minutes.
7. Although the BabyChecker will automatically turn off after about two minutes, slide the lb•oz/off/ kg switch to OFF. Then if the Start/Zero button is accidentally pressed, the batteries will not be used.

TIPS FOR USING THE BABYCHECKER

- When pressing the Start/Zero Button, do not touch any other part of the scale. Touching other parts might cause the scale to display negative figures such as -0.04 or -0.10 , therefore making the accuracy of the weight calculation invalid.
- When lifting or moving the scale, be sure to lift the entire scale from the bottom of the Scale Housing, not the Baby Basket.
- Do not put water on the scale or place the scale in a humid location.
- The scale is a precision instrument. Do not drop it, jump on it, or drop a heavy object onto its surface.
- Do not dismantle the scale.

CARE AND CLEANING INSTRUCTIONS

The scale is a precision instrument and specially adjusted and manufactured for accurate measuring. Use of the scale following these instructions is important.

1. Do not dismantle the Scale Housing.
2. The scale should not be operated in the presence of external vibration or direct sunshine.
3. Rapid change of ambient temperature may affect the accuracy of the scale for up to 2 hours.
4. Do not spill water or chemical liquid on or in the Scale Housing. Periodically, wipe the scale with a damp cloth. The Baby Basket should be cleaned using Cavicide® or another non-staining germicidal cleaner.
5. Keep the scale away from any heating appliances (oven, stoves, etc.).
6. Keep the unit in a cool, dustless and dry environment and remove batteries during extended storage.
7. When in storage, do not put anything on the Baby Basket.

SPECIFICATIONS

CAPACITY: 20 Kg or 40 lb.

DISPLAY ELEMENT: LCD.

POWER SOURCE: 6 Alkaline AA Batteries.

RATED VOLTAGE: 9VDC.

OPERATING TEMPERATURE: 15° – 35° C. or 59° – 95° F.

DIMENSIONS: Base – 10-1/2" length x 11-1/2" depth x 2-1/8" height.

Tray – 21-2/3" length x 10-3/4" depth x 3-1/8" height.

Base with Tray – 12-3/4" length x 11-1/2" depth x 5-1/4" height.

WEIGHT: with Tray – Approximately 6.5 lbs. or Approximately 3 Kg.

BATTERY LIFE: Approximately 100 hours.

SERVICE INSTRUCTIONS

■ Verification of Accuracy

Periodically, the scale should be tested to verify the accuracy of its weight measurement. This may be accomplished by using a 10 kg calibration weight, class F or better. These weights are available from Medela. Use the following instructions for verifying the scale's accuracy.

1. Switch the lb·oz / off / kg switch to kg.
2. Push the Start/Zero Button.
3. Within 60 seconds, gently place a 10 kg weight on the Baby Basket.
4. If the scale does not read between 9.96 kg – 10.06 kg, call Medela for return and repair information.


■ Battery Indicator

If Lo appears on the Display Panel, replace the batteries with new ones at once.

The image shows the 'Lo' indicator on a digital display panel. The characters 'L' and 'O' are rendered in a stylized, segmented font typical of LCD displays. The 'L' is on the left and the 'O' is on the right, both appearing as dark shapes against a light background.

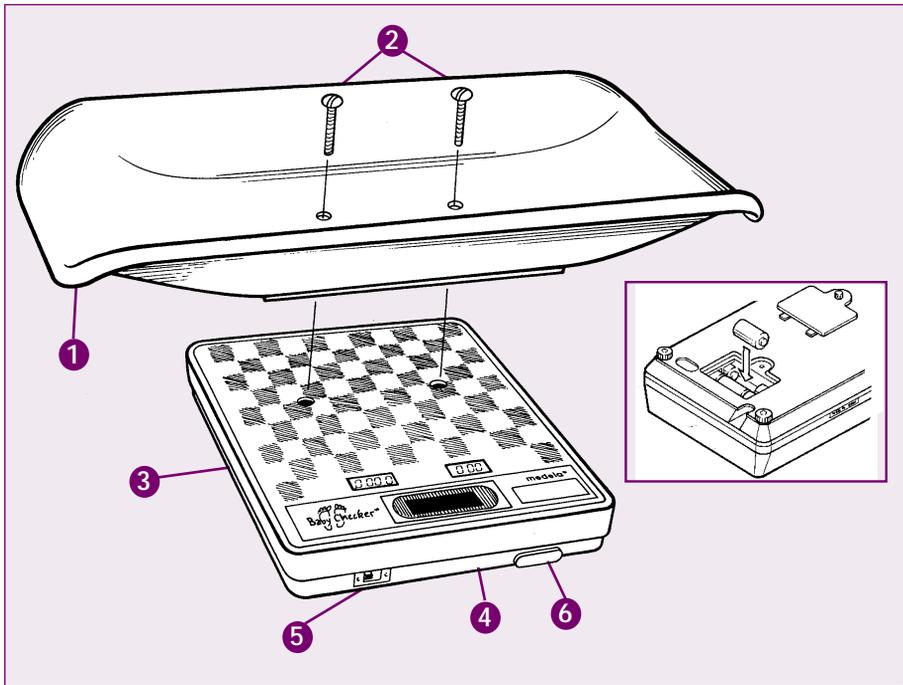If the batteries become too low, the Lo indicator will not display. If the scale does not turn on, replace the batteries to verify if this is the problem.

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©1997 Medela, Inc.
BabyChecker and natural solutions for beautiful beginnings... are trademarks of Medela, Inc.
Cavicide is a registered trademark of Micro-Aseptic Products, Inc.
Medela, Inc., P.O. Box 660, McHenry, IL 60050-0660 USA
Medela, Inc., P.O. Box/C.P. Box 131, Mississauga, Ontario, Canada, L4T 3B5
BabyChecker et des solutions naturelles pour débuter en beauté... sont des marques de commerce de Medela, Inc.